




Applied Scientific Instrumentation

Stage Inserts

283 x 110 mm Inserts

Inserts for ASI's MS-9400 XY Stage

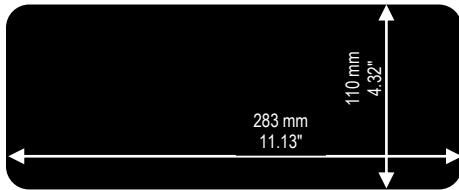

ASI's MS-9400 XY stage is designed for larger upright microscopes like the following:

- Leica: DMR series
- Nikon: FXA, SA, E800, E1000, 80i, 90i
- Olympus: AX70, BX50/60, BX51/61
- Zeiss: Axioskop, Axioskop II, Axioplan, AxioImager




Eight-Slide Insert for MS-9400 with Brass Sliders

I-9408-B

The I-9408-B slide insert accepts eight 25x75 mm (1" x 3") slides, and has individual spring-loaded brass slider restraints for easy loading and unloading. The unit is recessed to place the bottom of slides about 7.3 mm below the top of the insert. The slides can be easily preloaded into the insert before installing the insert into the stage. Multiple inserts can be used in this manner to screen a large number of samples.

Depth from Top of Insert: 7.3 mm
Overall Thickness: 8.0 mm

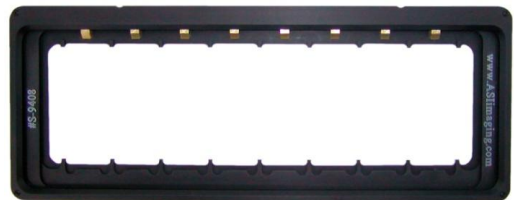


Eight-Slide Insert for MS-9400 with Spring Clips

I-9408-C

The I-9408-C slide insert accepts eight 25x75 mm (1" x 3") slides, and has individual spring clips for secure stability. The unit is recessed to place the bottom of slides about 7.3 mm below the top of the insert. The slides can be easily preloaded into the insert before installing the insert into the stage. Multiple inserts can be used in this manner to screen a large number of samples.

Depth from Top of Insert: 7.3 mm
Overall Thickness: 8.0 mm

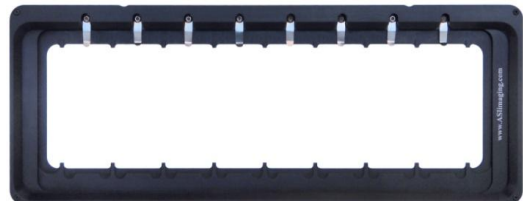

60 mm Three-Dish Insert for MS-9400

I-9460

The I-9460 insert accepts three Petri dishes with a maximum bottom outside diameter of 62.0 mm (2.44"). The unit places the bottom of the dish about 7.3 mm below the top of the insert. Surrounding cavity diameter is 67.1 mm (2.64") to accommodate a dish cover. Hole diameter is 57.0 mm (2.25"). Ease of access side slots. Spring clips for secure stability.

Depth from Top of Insert: 7.3 mm
Overall Thickness: 8.0 mm



We Create Solutions

Applied Scientific Instrumentation, Inc. ♦ 29391 W Enid Rd ♦ Eugene, OR USA 97402-9533
(541) 461-8181 ♦ (800) 706-2284 ♦ info@ASImaging.com ♦ www.ASImaging.com